


# CUNY test

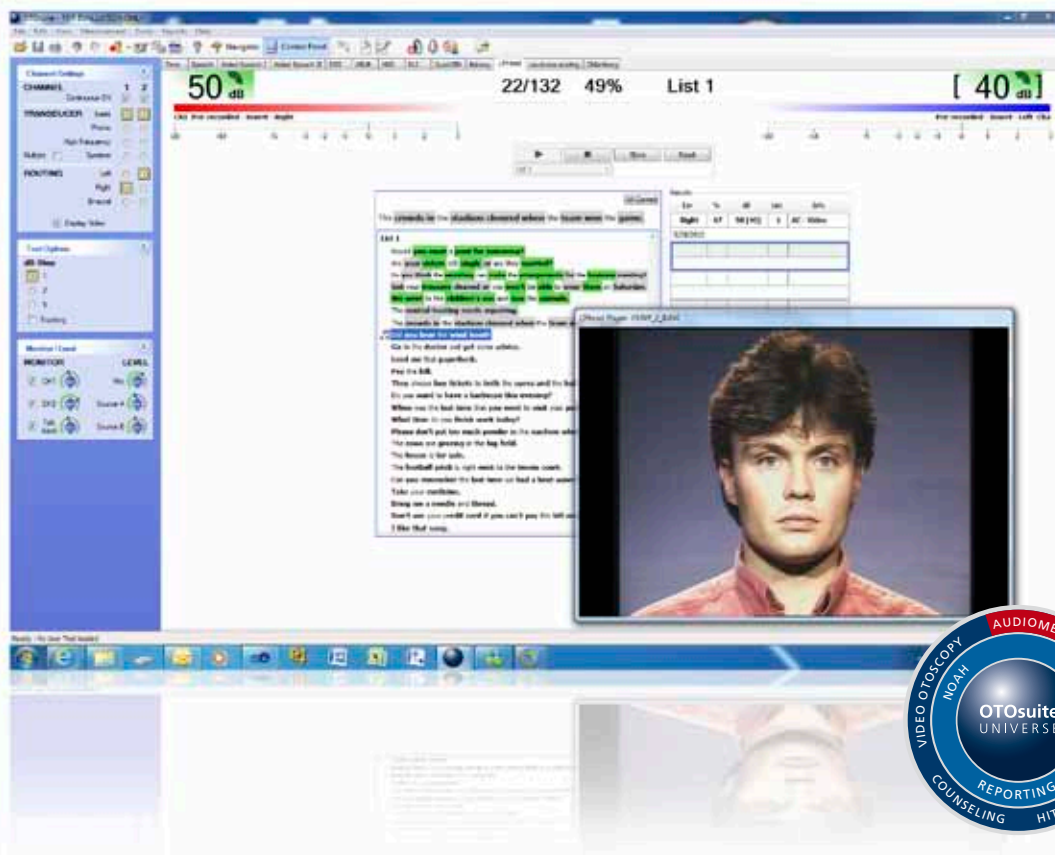

The original CUNY Sentence Test was developed by the City University of New York and consists of 72 lists with 12 sentences each. Each list contains 102 words and is scored for the total number of words correctly identified. The CUNY implementation in OTOsuite is a bimodal Lipread test recorded for use in the United Kingdom.

It is a test which is particularly used in cochlear implants the departments and used for assessing open-set sentence recognition. Before its integration into OTOsuite the words and sentences used for this test were recorded on tape and used by all cochlear implant centers. All of the tests are of a man's voice and played at the 70 Decibel range.

Using the CUNY test in OTOsuite gives you the following advantages:

- Easy navigation between the different tests performed (Tone and Speech audiometry and immitance)
- Reports with all the results
- A practical record of all the historic results and test comments in a single table
- A modern layout
- All sentence scoring details are stored and can be printed

## OTOsuite™

OTOSuite software is intuitive, easy to learn and easy to use. It supports the entire client journey, from diagnostics and counseling, to fitting and verification. Simple to use user test functionality ensures configurable, consistent and smooth workflows that can be customized to suit individual needs.

And you work in a single software environment across modules. The result is a better fitting experience – for you and for your client.

# Overview of audiometry tests in OTOsuite

	MADSEN Astera	AURICAL Aud
--	---------------	-------------



TEN	✓	✓
QuickSIN (US only)	✓	✓
CUNY (UK only)	✓	✗
Békésy	✓	✗
SAL	✓	✓
SISI	✓	✓
ABLB	✓	✓
Multi Frequency Weber	✓	✓
Oldenburger	✓	✓
Rinne	✓	✓
Tone Decay	✓	✓
DLI (Lüscher)	✓	✗
Stenger	✓	✓
MHA (Master Hearing Aid)	✓	✓
MLD (Masking Level Difference)	✓	✗

✗ Not Available    ✓ Available free of charge    ✓ Licensed Option

## Otometrics

Otometrics is the world's leading manufacturer of hearing and balance instrumentation and software. Over the last 50 years we have provided solutions ranging from newborn hearing screening applications and audiologic diagnostics to comprehensive hearing instrument fitting and balance testing.

However, there is more about Otometrics than just instrumentation. The value of Otometrics as a partner is about our commitment to our users and their clients. The reason why major sites around the world have chosen us as their preferred partner is more about trust in basing their mission-critical activities on what we provide. Our strong worldwide team of experts is determined to carry this legacy into the future.

 [facebook.com/otometrics](https://facebook.com/otometrics)  
 [twitter.com/otometrics](https://twitter.com/otometrics)

**GN Otometrics, Europe.** +45 45 75 55 55. [info@gnotometrics.dk](mailto:info@gnotometrics.dk)  
**GN Otometrics, North America.** 1-800-289-2150. [sales@gnotometrics.com](mailto:sales@gnotometrics.com)  
[www.otometrics.com](http://www.otometrics.com)



**otometrics**  
MADSEN · AURICAL · ICS